


SaniSure CLAMP

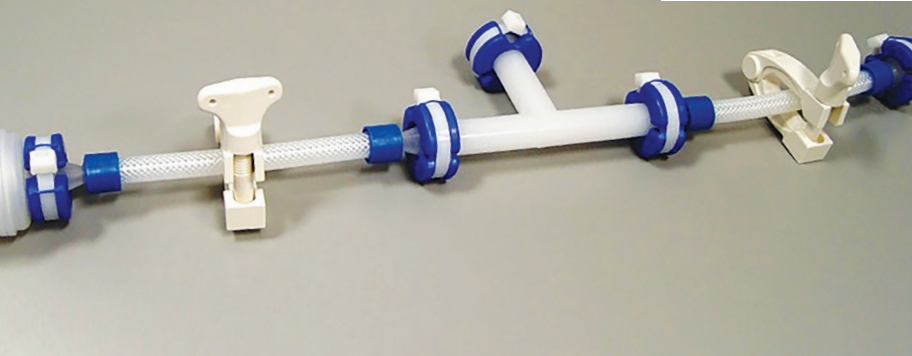

The SaniSure CLAMP is the very first sanitary clamp designed to offer a fast, safe and secure way to connect plastic to plastic, plastic to stainless steel and stainless steel to stainless steel for sanitary fittings. The SaniSure CLAMP works by providing conforming compression on tri-clamp fittings, which results in an even clamping force on the entire fitting surface — something a traditional “rigid” clamp design can’t achieve.

- Not Operator Dependent
- Validatable
- Improve Safety
- Reduce Setup Time
- Meets USP Class VI
- Patent Pending Design

FEATURES

BENEFITS

360-Degree Compression

Prevents leakage caused by over-tightening of traditional clamps; Reduces gasket intrusion caused by traditional two-point load clamps; Ability to handle a wide range of pressures up to 100 psi

Low-Profile and Lightweight

No sharp edges, reducing risk of penetrating bag assemblies; Quick and easy to secure for operators; Reduces shipping weight costs

Cable Tie Locking Mechanism

Reduces risk associated with carpal tunnel and wrist injuries; Precision calibrated tension ensures consistent and repeatable results; Tamper proof and safe

Controlled Compression Clamp


Connect plastic to plastic fittings; Connect plastic to BPE stainless steel fittings; Connect BPE stainless steel to BPE stainless steel fittings

Single Use or Multiple Use

Cost effective; Eliminates cleaning validation; Reduces installation time; Durable enough to withstand repetitive use

Installation Procedures

Simple three-step installation

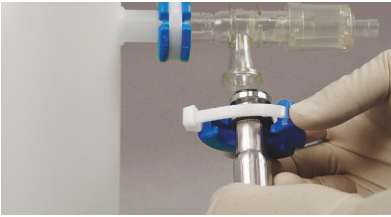

Position SaniSure Clamp around tri-clamp fitting



Insert cable tie end and slide clamp closed



Finish closure by using SaniSure installation gun to calibrated tension



Controlled Compression

SaniSure clamps work by providing conforming compression on tri-clamp fittings. The result is even clamping force on the entire fitting surface that no traditional “rigid” clamp design can achieve. The clamps work equally well on disposable and BPE stainless tri-clamp fittings in any combination.

Applications

- Bioreactors
- Single-Use Assemblies and Rigging Sets
- Plastic Filter Housings and Capsules
- Shipping Containers
- Contract Manufacturing
- Lock-Out, Tag-Out Critical Components
- Skid Systems
- Chromatography Columns
- Process Containers

Sterilization

Autoclave

Gamma Irradiated to Sterility Assurance Level 10⁻⁶

Part Number	Description	WP psi at 70°F	Material of Construction	Tri-Clamp Size
SS-NYB-MINI	3/4" Tri-Clamp	150	Nylon 6/6	Mini
SS-NYB-150	1 & 1 1/2" Tri-Clamp	100	Nylon 6/6	1 1/2"
SS-MK-9	Installation Gun			All Sizes
SS-MK-9-CAL	Calibrated Installation Gun			All Sizes
SS-CUT-170	Cable Tie Removal Tool			All Sizes
SS-CAB-TIE-750	Cable Tie for 3/4" Clamp		Nylon 6/6	Mini
SS-CAB-TIE-1500	Cable Tie for 1 1/2" Clamp		Nylon 6/6	1 1/2"



HENRICUSKADE 119A | 2497 NB THE HAGUE | THE NETHERLANDS | T +31 (0)85 043 31 10
F +31 (0)84 714 05 49 | E INFO@ROMYNOX.NL | I WWW.ROMYNOX.NL

IMPORTANT: It is the user's responsibility to ensure the suitability and safety of SaniSure, Inc. products for all intended uses.

For a period of 90 days from the date of delivery, SaniSure, Inc. warrants this product to be free from defects in materials and workmanship. Our only obligation will be to replace any portion proving defective, or at our option, refund the purchase price thereof. User assumes all other risk, if any, including the risk of injury, loss or damage, direct or consequential, arising out of the use, misuse or inability to use this product. THIS LIMITED WARRANTY IS IN LIEU OF THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR PARTICULAR PURPOSE AND ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED. No deviation is authorized.

SaniSure, Inc. assumes no obligations or liability for any advice furnished by it or for results obtained with respect to those products. All such advice is given and accepted at the user's risk. SaniSure® is a registered trademark of SaniSure, Inc. All rights reserved.